

ST216

Electrosonic Care+

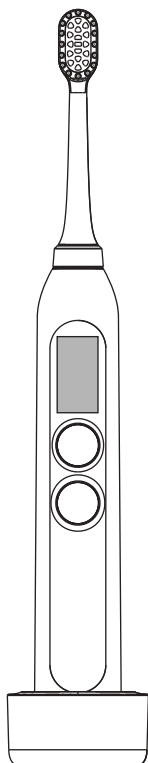
EN	Electric sonic toothbrush - user manual	4
PL	Szczoteczka elektryczno-soniczna - instrukcja obsługi	8
DE	Elektrische Schallzahnbürste - Benutzerhandbuch	12
BG	Електрическа звукова четка за зъби - ръководство за потребителя	16
CS	Elektrická sonická zubní kartáčka - uživatelská příručka	20
DA	Elektrisk sonisk tandbørste - brugervejledning	24
EL	Ηλεκτρική οδοντόβουρτσα - εγχειρίδιο χρήσης	28
ES	Cepillo de dientes eléctrico sónico - manual de usuario	32
ET	Elektriline sonic-harja - kasutusjuhend	36
FI	Sähköinen sonic-hammasharja - käyttöohje	40
FR	Brosse à dents électrique sonique - manuel d'utilisation	44
HU	Elektromos szonikus fogkefe - felhasználói kézikönyv	48
IT	Spazzolino elettrico sonico - manuale d'uso	52
LT	Elektrinis soninis dantų šepetėlis - naudotojo vadovas	56
LV	Elektroniska soniska zobu suka - lietotāja rokasgrāmata	60
NL	Elektrische sonische tandenborstel - gebruikershandleiding	64
NO	Elektrisk sonisk tannbørste - brukerveiledning	68
PT	Escova de dentes elétrica sónica - manual do utilizador	72
RO	Periuță de dinți electrică sonică - manual de utilizare	76
SK	Elektrická sonická zubná kefka - používateľská príručka	80
SV	Elektrisk sonisk tandborste - bruksanvisning	84


1

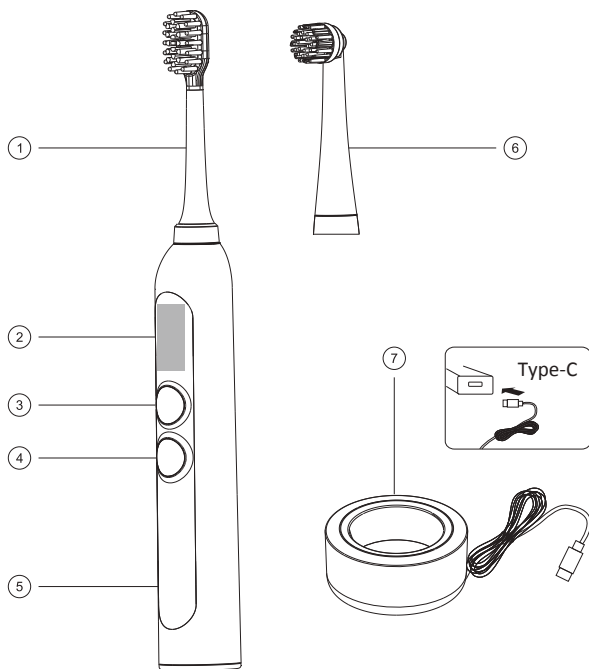
2

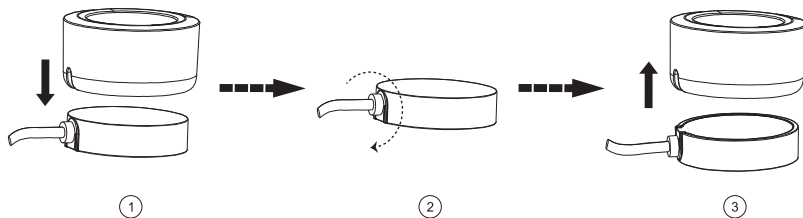
3

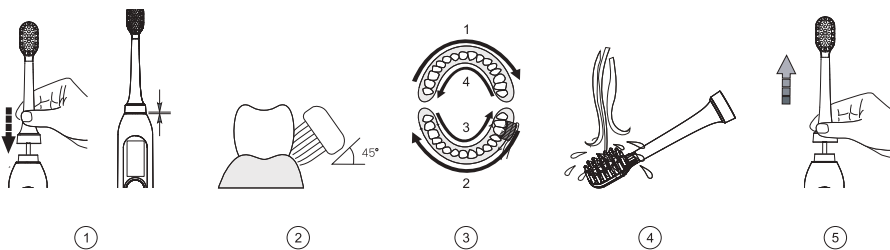
4

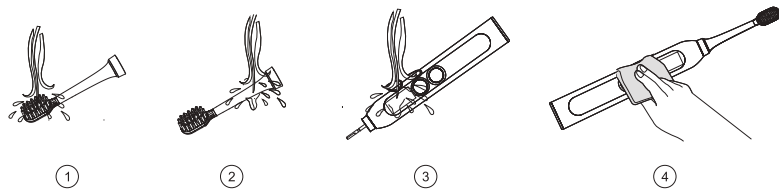5

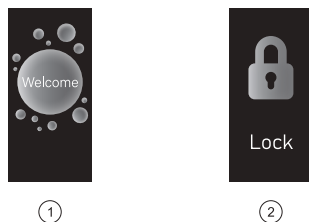
6

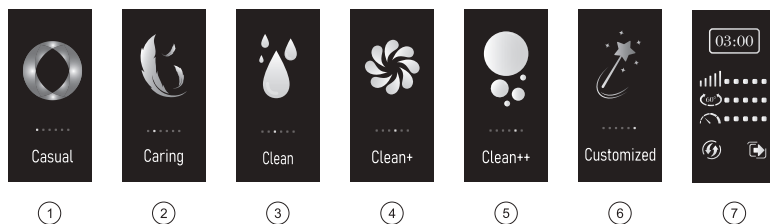7

8

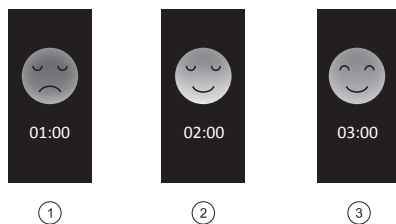Safety instructions & notice

Please read this manual carefully to understand better how to use the product. Some essential safety tips, including the following, shall be observed during use.

- Do not expose the device or batteries to excessive temperatures.
- Please note that battery-powered electrical devices are at risk of short-circuiting.
- If you notice a battery leak, stop using the device immediately and dispose of it in accordance with local environmental regulations and guidelines. Do not dispose of a leaking battery with normal household waste.
- Do not place the charging parts of this product in a damp environment or in water.
- Do not charge the device if there is liquid on the charging interface.
- To ensure the proper functioning of this product, use it away from magnetic objects or in an environment with a strong magnetic field.
- This device contains a battery that cannot be replaced. When the battery is depleted, dispose of the device in accordance with environmental recycling instructions.
- To ensure effective battery charging, charge the battery within a temperature range of 0°C to 45°C to avoid incomplete or poor charging.
- Check the device regularly for damage. If damage is found, stop using it immediately and replace it in a timely manner.
- Do not hit or drop the product. If the product is damaged, especially the brush handle and brush head, do not continue to use it.
- Do not modify or repair this product yourself.
- Do not use chemicals such as paint thinner to clean this product.
- Do not use the product in water deeper than 1 meter for more than 30 minutes.
- When using the product for the first time, you may experience temporary discomfort in your teeth or slight bleeding from your gums. This is because your gums are not yet accustomed to intensive cleaning. We recommend adjusting the product to your preferred gentle setting. These symptoms usually disappear within 1-2 weeks. If they persist for longer, discontinue use immediately and consult your dentist.
- This product is intended for cleaning teeth, gums and tongue only. Do not use it for any other purpose.
- Users of dentures, orthodontic appliances, dental implants, etc. should use the product under the supervision of a dentist; in case of significant bleeding, discomfort, pain or other health problems, consult your dentist or doctor.
- The product is not suitable for children, physically disabled persons, persons with certain sensory and mental disorders, and persons with a serious lack of experience. Use this product safely under supervision and understand the potential risks.
- Do not share the brush head with others.
- It is recommended to replace the brush head every three months. If the brush head is damaged or deformed, replace it immediately with an officially certified replacement part.
- This product may only be used with an electric toothbrush from this company and a special brush head. We are not responsible for any consequences resulting from the use of unofficial, uncertified products with the brush head.

Technical parameters

Product name: Vibrosonic / Rotation electric toothbrush with screen

Model: ST216

Charging input: 5 V - 1 A

Rated power: 5 W

Battery capacity: 2000 mAh

Charging time: Approximately 16 hours

Waterproof rating: IPX7

Main unit description

Illustration no. 1

1. Vibrosonic brush head
2. LCD screen
3. On/Off button
4. Modes button
5. Handle
6. Rotating brush head
7. Charging station

Battery level

1. Battery charge indicator: (0%–100%).
2. The remaining charge icon will be displayed each time the toothbrush is turned on or the brush head is reset.

3. During charging, the current charge level icon is displayed, and when fully charged, the 100% charge icon is displayed.
4. During discharge protection, the motor stops working. The battery icon does not show power, and the toothbrush turns off immediately.
5. If the toothbrush cannot be charged after placing it on the charging base, the charging plate may be installed upside down. Turn the plate upside down before charging. (illustration no. 2)

Operation mode

1. Automatic brush head reset: When the handle is lifted and touch is detected, the brush head automatically resets to the centre position and the display enters standby mode.
2. Welcome screen and locked screen. (illustration no. 5)
3. When the device is turned on using the power button or when the brush head is reset by lifting it, the screen displays a welcome message, the current battery level and the notification 'End of playback'. The system will display the current operating mode and enter standby mode. Pressing the power button again will activate the motor. Pressing the power button again will reset the motor. When the vibration stops, the current score is displayed on the screen. When the screen is turned off, the toothbrush enters the off state. In standby mode, after 5 seconds of inactivity, the toothbrush automatically enters the off state.
4. Lock and unlock functions: In off, standby or active mode, press and hold the mode button for more than 5 seconds to enter lock mode. Once successfully locked, the lock icon will remain lit for 5 seconds. Similarly, in locked mode, press and hold the mode button for more than 5 seconds to unlock. During unlocking, the lock icon will flash intermittently; when complete, it will display the current operating mode and enter standby mode.

Vibration mode

1. When the brush head is attached to the handle or no brush head is in use, the default mode is vibration.
2. Vibration mode interface. (illustration no. 6, points 1–6)
3. By pressing the mode switch button, you can switch between six modes: Casual, Caring, Clean, Clean+, Clean++, and Customized.
4. Mode customization interface. (illustration no. 6, point 7)
5. In standby/working mode, press and hold the on/off button for 5 seconds to enter the customized settings interface; after 8 seconds of inactivity, the customized settings interface will automatically close; select the return function to exit the settings interface.
6. Definitions of icons in the custom settings interface: intensity, angle, speed, reset, return. After selecting a function, the icon or number will flash. Each press of the on/off button adds or confirms a function.
7. Brushing time is optional: 1 minute, 1 minute 30 seconds, 2 minutes, 2 minutes 30 seconds, 3 minutes.
8. The reset function restores the factory default settings with a force of 5, an angle of 1 and a speed of 1. Each addition of a level means the current intensity, speed and angle.

Rotation mode

1. When the rotating brush head is put on the handle, it automatically switches to rotation mode.
2. Rotation mode interface. (Illustration no. 7)
3. There are six rotation modes, which can be switched by pressing the gear button in a loop: Soft spin, Auto spin, Ultra spin, Soft osc, Auto osc, Ultra osc.

Display rating

1. At the end of each brushing cycle, the current brushing effect is scored. Before shutdown, the current brushing cycle score is displayed on the screen to remind users whether the brushing time is sufficient. (Illustration no. 8)

Automatic functions

1. Automatic timer function: In each mode, the brushing time is set to 2 minutes by default. The custom mode is set to 2 minutes by default. During the timer cycle of the selected mode, changing the mode will not reset the timer (the time is cumulative), and the toothbrush will continue to operate according to the timer cycle of the original mode.
2. Interval reminder function: In each mode, the toothbrush reminds the user to change the brushing area every 30 seconds.
3. Auto memory function: After using any mode for more than 5 seconds, the toothbrush will remember the selected mode. The next time the power is turned on, it will automatically start operating in the remembered mode.

Daily use

Illustration no. 3

1. Attach the brush head to the handle as shown in the diagram, leaving approximately 0.6 mm between the brush head and the end of the handle.

2. Rinse the brush head and apply an appropriate amount of toothpaste. Position the bristles at a slight angle to the gum line and then place them on your teeth.
3. Gently press down and slowly move the brush head, making small movements between your teeth. Avoid using excessive force.
4. To clean your entire mouth more effectively, divide it into four quadrants. Start with the first quadrant (outer upper teeth) and brush each quadrant for 30 seconds. Every 30 seconds, brushing will pause to remind you to move on to the next quadrant. The two-minute smart timer will automatically turn off when the brushing time is up.
5. After brushing all four quadrants, you can turn the device back on to brush the chewing surfaces of your teeth and areas with tartar build-up.
6. Clean the brush head and handle by rinsing them with water. Gently tap the brush head to ensure that no food debris remains in the bristles.
7. Remove the brush head in the direction indicated in the diagram and shake off any excess water.

Toothbrush cleaning

Illustration no. 4

1. Gently remove the brush head from the handle.
2. Rinse the brush head and bristles with clean water to remove food debris and toothpaste foam.
3. It is recommended to remove the brush head from the handle and rinse the connection point at least once a week to prevent the growth of bacteria.
4. Wipe the surface of the handle with a damp cloth to remove dirt and debris.
5. Use a dry soft cloth to dry the handle.

Troubleshooting

If the toothbrush malfunctions, such as the brush head failing to return to its original position, rendering the button functions inoperable, please immediately stop using it and contact customer service for after-sales handling. We will provide you with professional solutions and services promptly.

Warranty terms and conditions

1. The product comes with a 24-month manufacturer's warranty, carried out by sole distributor: Morele.net sp. z o.o., al. Jana Pawła II 43b, 31-864 Kraków, NIP: PL9451972201.
2. Complaints / warranty repairs should be reported through the warranty / service form available on the morele.net website.
3. The warranty covers only latent defects attributable to the manufacturer.
4. The claim shall be processed by the point of sale or technical service under the warranty if the Customer delivers
 - a. a valid proof of purchase with the date of sale,
 - b. the reported product or defective part. In the case of a remote purchase, the warranty card is valid solely with proof of purchase (receipt/invoice).
5. Warranty exclusions:
 - a. mechanical damage and related defects,
 - b. damage and defects resulting from misuse and inappropriate storage, assembly, and maintenance,
 - c. damage and wear of such consumable elements,
 - d. As per the manual, activities related to assembly or maintenance are the user's obligation.
6. The warranty shall not apply when:
 - a. it is past its validity period,
 - b. the Customer had repaired or modified the product using non-genuine parts,
 - c. the defect resulted from incorrect installation or a failure to follow the rules of correct operation by the manual,
 - d. the product is used for purposes other than household use,
 - e. damage occurred during transport.
7. The Customer is entitled to the following compensation free of charge:
 - a. product repair,
 - b. product replacement,
 - c. reduced price,
 - d. termination of a contract and full cost compensation.
8. To report a claim
 - a. Present the product or its reported part.
 - b. Present the proof of purchase with the seller's name and address, date and place of purchase, and type of product, or a valid warranty card with a stamp of the point of sale.
 - c. If the product is delivered dirty, the technical service may refuse to accept it or clean it at the Customer's cost after receiving their written approval.

9. If the claim is accepted, the device is repaired or replaced with a new one, or the Customer is reimbursed. Transport costs from the manufacturer's service cover service.
10. We are not liable for damage caused by improper use of the device. Costs of ingredients wasted as a result of improper operation of the device, lost profits, etc. are not refundable.

Utilization

Handling of end-of-life waste from electrical and electronic equipment (WEEE legislation).

The marking on the packaging, accessories, or instruction manual indicates that end-of-life equipment and attached electronic accessories should not be placed in municipal waste containers with other household waste.

Improper handling of electronic waste may result in contaminants entering the environment or health risks due to inappropriate disposal. Separate this waste from other waste and recycle it so that the raw materials can be reused. For recycling, they are allowing the reuse of raw materials.

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